

18.5 cu Large Capacity Combination Blanket / Fluid Warmer

Item# LB0581CEEMS

Warm your blankets and fluids for busy emergency rooms, operating rooms and labor and deliver units in the convenient, reliable and affordable combination warmer. Separate controls ensure that blankets and fluids are safely and accurately heated to recommended temperatures. Improve recovery time by efficiently warming patients from admission to discharge, reducing surgical site infections.

Features:

- The durable, heavy-duty stainless steel blanket warmers quickly and efficiently heat blankets and fluids to precise temperatures using the digital control, ensuring your warmed fluids meet healthcare guidelines for IV and irrigation fluids.
- Ensure each bag is safely warmed to recommended temperatures through the WarmSafe zone heating system. WarmRight technology provides the right temperature for the right product to meet healthcare guidelines for IV and irrigation fluids.
- Dual-chamber warming appliance
- Stainless steel exterior and interior
- Doors are fully gasketed with a magnetic closure
- WarmSafe Zone Heating: Heat is evenly balanced throughout the chamber to provide safe fluid warming.
- WarmWatch - IV bag Temperature Reader: Know the temperature of every bag of fluid before administering to patients.
- Convenient: Easily stock and access blankets and fluids from the front and back with the optional pass-through door.
- Durable: The affordable, rugged stainless steel cabinets are made to last and work just as hard as you do.
- Mobile: Rugged casters are designed to easily transport warmed fluids and blankets throughout the facility.



Toll Free 1.877.706.4480 | Local/ International 713.706.4480 | Fax 713.706.4161
www.MSEC.co

- Glass Window: The glass window allows staff to take an immediate inventory of the amount of warmed fluids and blankets available
- Service: Our worldwide support team is dedicated to providing quality equipment. We walk with you every step of the way to maximize your return on investment.

Control:

- Blanket chamber has an adjustable temperature range of 32°C - 71°C (90°F - 160°F)
- Fluid chamber can be programmed to warm either irrigation fluids (IRR) or injection fluids (INJ), with separate temperature ranges provided depending on the choice selected:
- IRR temperature range: 32°C to 66°C (90°F to 150°F)
- INJ temperature range: 32°C to 40°C (90°F to 104°F)
- Accurate temperature
- IRR +0/-1.67°C (-3°F)
- INJ +0/-1.10°C (-2°F)
- Operate in Celsius or Fahrenheit
- LED display
- Built-in speaker for audible feedback
- Integrated controller lock-out

Specifications:

- Capacity: 30-40 blankets; 30 1-liter bottles or 24 1-liter bags
- Interior: 18.5 ft³
- Size: 74.8 in H; 30 in W; 28 in D



40

30-40 blankets

30

30 (1) liter bottles/24 (1) liter bags



Adjustable temperature range of

32°-71°C (90°-160°F) adjustable blanket mode

32°-40°C (90°-104°F) adjustable injection mode

32°-66°C (90°-150°F) adjustable irrigation mode

Accurate temperature


IRR +0/-1.67°C (-3°F)

INJ +0/-1.10°C (-2°F)



Additional features:

- Safety: The controller stores all settings and resumes operation using these settings in the event of a power failure. When the power is restored, the controller sounds an alarm. The on/off status indicator flashes.
- Safety: A warming shut-off system, separate from the controller, prevents overheating.
- Convenience: Access point at rear panel allows the addition of data logging or temperature management hardware
- Manufacturer warranty: 1 year labor & 3 year parts

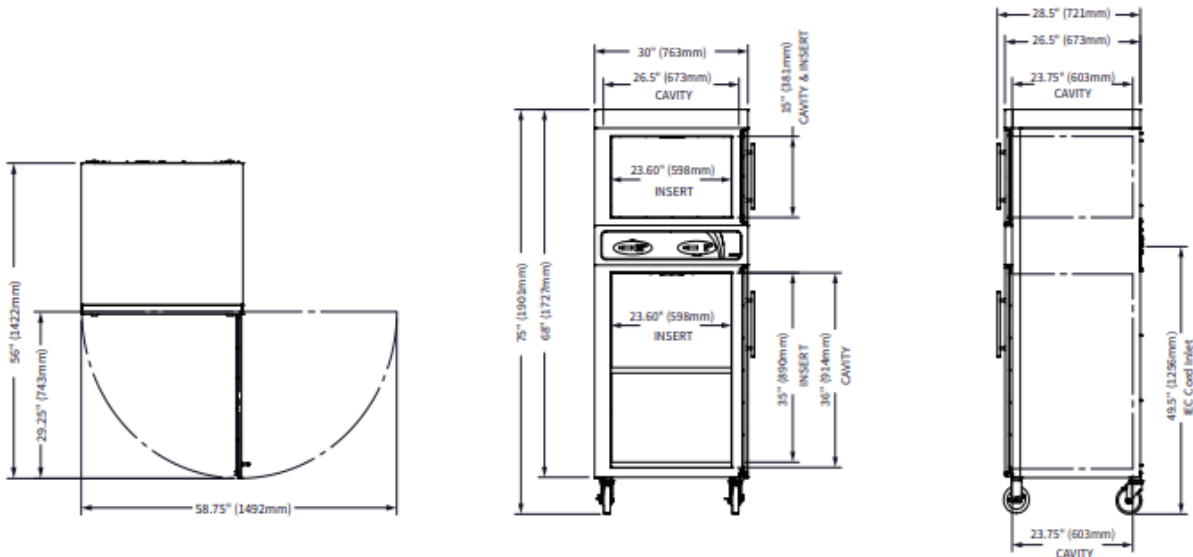

Top: 1"
Left: 2"
Right: 2"
Back: 4"

CLEARANCE



Exterior (H x W x D)

75" (1905mm) x 30" (762mm) x 28.5" (724mm)

Interior (H x W x D)

51" (1295mm) x 26.5" (673mm) x 23.75" (603mm)

Net Weight

457lbs (207kg)

Ship Dimensions (L x W x H)*

38" x 35" x 82"

Ship Weight*

528 lbs

*Domestic ground shipping information. Contact factory for export weight and dimensions.



ELECTRIC

V	Ph	Hz	A	kW	Plug Configuration
120	1	50/60	11	1.3	NEMA 5-20P; 20A 125V Plug; Hospital Grade
220	1	50/60	5.5	1.2	CEE 7/7 220-230V Plug
230	1	50/60	5.6	1.3	BS 1363 Plug; CEE 7/7 220-230V Plug

*Other international plugs are available, contact the manufacturer for more information.

**Electrical connections must meet all applicable federal, state and local codes.



Toll Free 1.877.706.4480 | Local/ International 713.706.4480 | Fax 713.706.4161
www.MSEC.co