

Fluid Warming



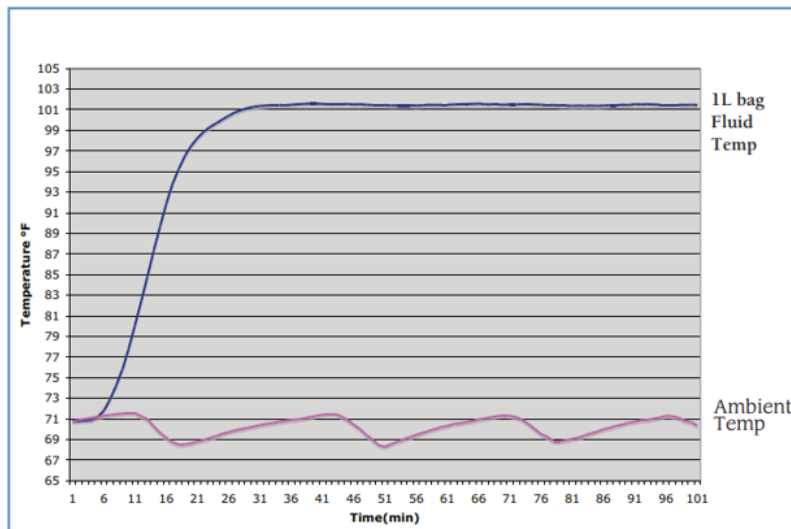
- Modular configuration Wall mount: 1-6 bag capacity Countertop: 1-3 bag capacity
- Pole mount: 1-3 bag capacity
- Accommodates 0.5L, 1L, 2L & 3L bags
- Anodized aluminum heat plate
- Lightweight, medical grade plastic housing
- Intelligent bag sensors
- Dual temperature sensors
- Heats to 40°C (104°F) in 30 minutes or less*
- *3L bags require slightly longer heating time
- Accuracy of +0°/-2°C (+0°/-3°F)
- Seismic Horizontal Wall Mount Bracket Kit is OSHPD approved. OPA-0698-10
- 3L Bag Wall Mounted Tilt Kit is OSHPD approved OPA-0060-13



Toll Free 1.877.706.4480 | Local/ International 713.706.4480 | Fax 713.706.4161
www.MSEC.co

Modular Fluid Warmer

The fluid warmer quickly warms and maintains the temperature of injection/intravenous and irrigation solutions prior to their use. The specially contoured warming module cradles solution bags in 0.5-, 1-, 2- & 3-liter sizes. Three-liter bag capacity available with an additional adapter. Individual units cannot be stacked in the field. Multiple cavity units are available from the factory in the configurations below. The unit is controlled by one (1) power switch and individual electronic controls with L.E.D. display for each cavity. The control can easily be set to display temperatures in Celsius or Fahrenheit. A sensor in the heating plate detects the presence of a bag and engages the heating mechanism to quickly begin warming the fluid. Two (2) temperature sensors work in unison to precisely and continuously read the temperature of the bag and another sensor monitors the plate temperature. A green ready light will illuminate when fluid is within $\pm 2^{\circ}\text{C}$ ($\pm 3^{\circ}\text{F}$) of the set point temperature. The heater will reengage as necessary to maintain the temperature within $\pm 2^{\circ}\text{C}$ ($\pm 3^{\circ}\text{F}$) of the set point. The electronic control monitors the length of time the bag has been held at temperature, beginning when the bag reaches set point temperature. A status button will display the time the fluid has been held at temperature.



Available Configurations:

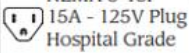

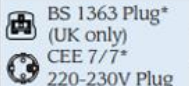

NOTE: Individual warmers cannot be stacked in the field. Multiple cavity warmers are available from the manufacturer in the configurations listed below. 3-liter bag capacity is available with an additional adapter. ivNow-1 (one (1) bag capacity) can be placed directly on a countertop, mounted on a mobile equipment pole stand, mounted on a wall using brackets, or mounted to a 3-liter bag tilt kit*. Two (2) bag capacity can be placed directly on a countertop, mounted on a wall using mounting brackets, or mounted to a 3-liter bag tilt kit*.


Three (3) bag capacity must be mounted on a heavy-duty mobile equipment pole stand, on a wall using mounting brackets, or mounted to a 3-liter bag tilt kit*. Four (4) bag capacity must be mounted on a wall using mounting brackets*. Five (5) bag capacity must be mounted on a wall using mounting brackets*. Six (6) bag capacity must be mounted on a wall using mounting brackets*. *Mounting hardware sold separately.



Toll Free 1.877.706.4480 | Local/ International 713.706.4480 | Fax 713.706.4161
www.MSEC.co


Modular Fluid Warmer

Specs	ivNow-1	ivNow-2	ivNow-3	ivNow-4	ivNow-5	ivNow-6	Plug	Certification
Height (wall unit)	8.0" (203.1mm)	13.4" (340.5mm)	18.9" (478.8mm)	24.2" (615.3mm)	29.6" (752.7mm)	35" (890.1mm)	—	—
Width (wall unit)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	—	—
Depth (wall unit)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	—	—
Cavities	1	2	3	4	5	6	—	—
Net Weight	8 lbs (4 kg)	14 lbs (6 kg)	20 lbs (9 kg)	25 lbs (11 kg)	31 lbs (14 kg)	37 lbs (17 kg)	—	—
Ship Weight	11 lbs (5 kg)	16 lbs (7 kg)	23 lbs (10 kg)	28 lbs (13 kg)	34 lbs (16 kg)	40 lbs (18 kg)	—	—
120V							 NEMA 5-15P 15A - 125V Plug Hospital Grade	
AMPS//kW	1.3//0.15	2.5//0.3	3.8//0.45	5.0//0.6	6.3//0.75	7.5//0.9		
Hz//PH	60//1	60//1	60//1	60//1	60//1	60//1		
230V							 BS 1363 Plug* (UK only) CEE 7/7* 220-230V Plug	
AMPS//kW	0.7//0.15	1.3//0.3	2.0//0.45	2.6//0.6	3.2//0.75	3.9//0.9		
Hz//PH	50/60//1	50/60//1	50/60//1	50/60//1	50/60//1	50/60//1		

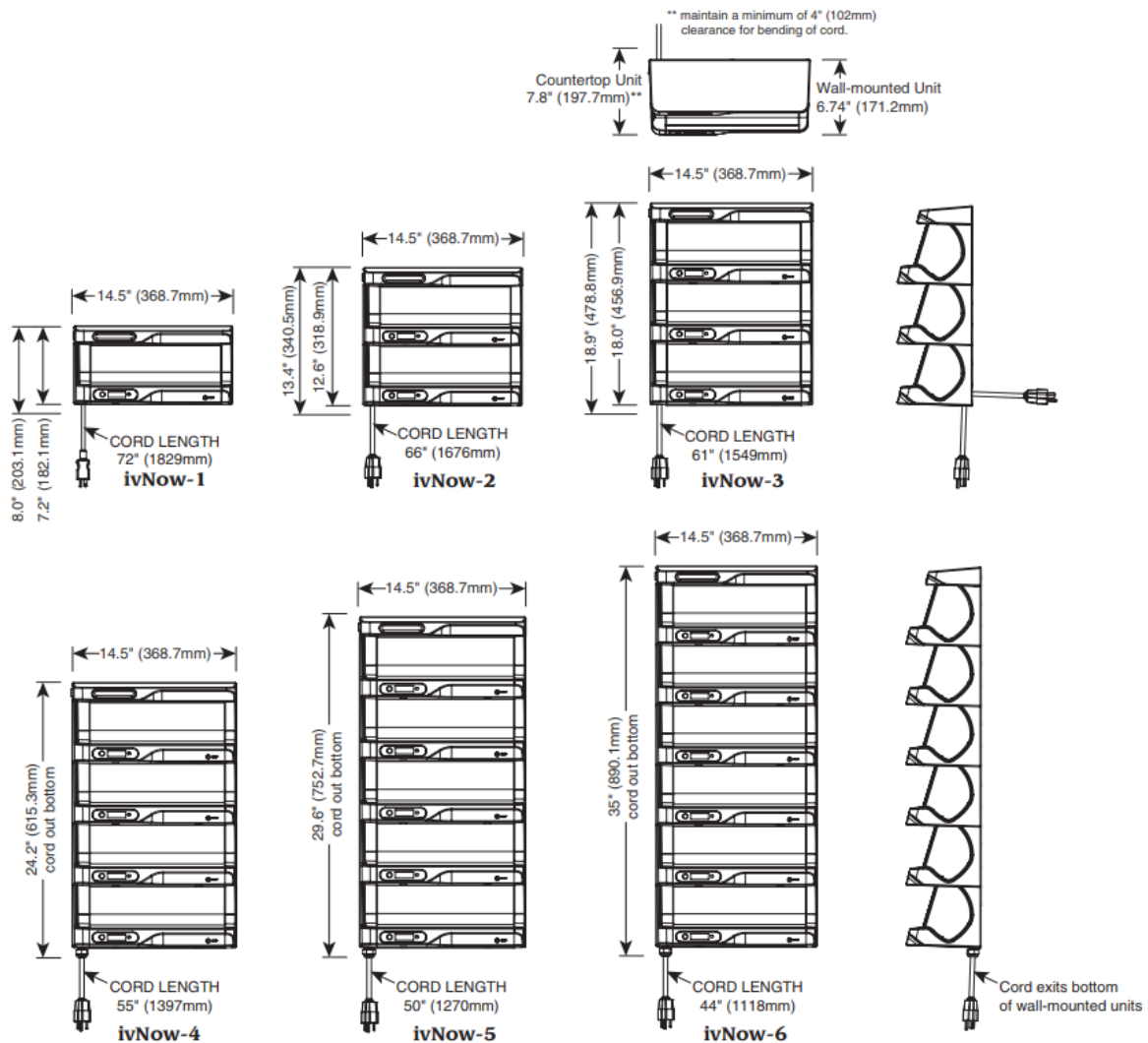


ACCESSORIES		ivNow-1	ivNow-2	ivNow-3	ivNow-4, -5, -6
3L Bag Tilt Kit	COUNTERTOP	5014126	5014127	5014128*	—
	WALL MOUNTED	5014129	5014130	5014131*	—
	WALL MOUNTED OPA-0060-13	5017700	5017701	5017702*	—
GCX Light Duty Roll Stand for Devices & Bracket Kit		mount to GCX equipment stand - stand included (cord exits bottom of unit) 5012246	—	—	—
GCX Heavy Duty Roll Stand for Devices & Bracket Kit		—	mount to GCX equipment stand - stand included 5013312	mount to GCX equipment stand - stand included 5013312	—
Bracket Kit for Equipment Stand		mount to equipment stand or mast - stand not included (cord exits bottom of unit) 5012735	—	—	—
Bracket Kit for Horizontal Wall Mount		two (2) GCX horizontal rails included 5012242	two (2) GCX horizontal rails included 5012242	two (2) GCX horizontal rails included 5012242	(cord exits bottom of unit) two (2) GCX horizontal rails included 5012242
Bracket Kit for Seismic Horizontal Wall Mount OPA-0698-10		one (1) GCX horizontal rail included 5016601	one (1) GCX horizontal rail included 5016601	two (2) GCX horizontal rails included 5016602*	two (2) GCX horizontal rails included 5016602*
Vertical Wall Channel Mount & Bracket Kit		one (1) GCX vertical rail included 5012241*	one (1) GCX vertical rail included 5012241*	one (1) GCX vertical rail included 5012241*	one (1) GCX vertical rail included 5012241*
Bracket Kit for Vertical Wall Channel Mount		GCX vertical rail not included 5012288*	GCX vertical rail not included 5012288*	GCX vertical rail not included 5012288*	GCX vertical rail not included 5012288*
Bracket for Mounting to Harloff Anesthesia Cart		1015253	1015254	1015255	—

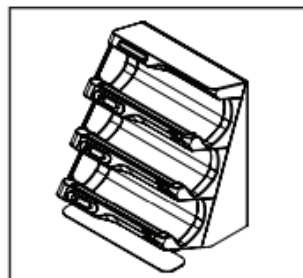
* See equipment diagram below



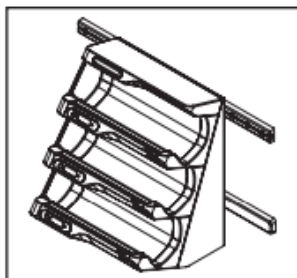
3L Bag Tilt Kit for countertop units (5014128 shown)



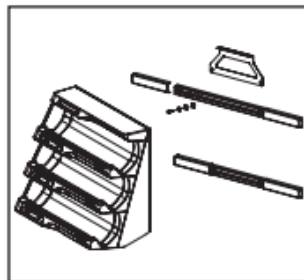
3L Bag Tilt Kit for wall mounted units (5014131 shown)



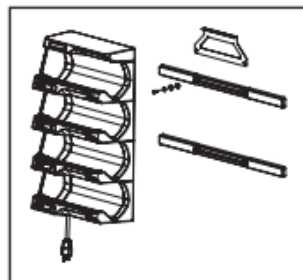
3L Bag Tilt Kit for wall mounted units (5017702 shown)
OPA-0060-13



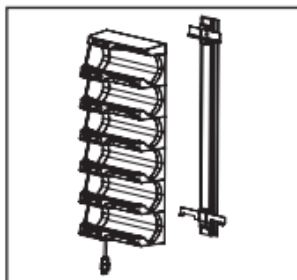Bracket Kit for Seismic Horizontal Wall Mount 5016602 (ivNow-3 to -6)
OPA-0698-10



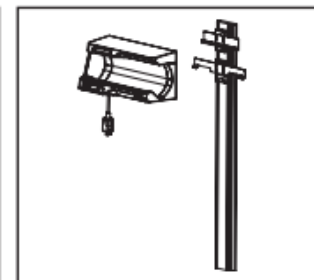(ivNow-6 shown above)
Vertical Wall Channel Mount Bracket Kit 5012241
Brackets only 5012288



(ivNow-1 shown above)



Toll Free 1.877.706.4480 | Local/ International 713.706.4480 | Fax 713.706.4161
www.MSEC.co